

Installation of Apache, PHP, Zend optimizer and MySQL in Windows


- I. Chemoventory
- II. Apache
- III. PHP
- IV. Zend
- V. MySQL

Please download latest stable versions of Apache 2.2, PHP, MySQL, Zend and MySQL.

I. Copy chemoventory in folder c:/chemoventory/

These are the folders you should see and about 130 files.

You can choose different directory, but make appropriate changes when required.

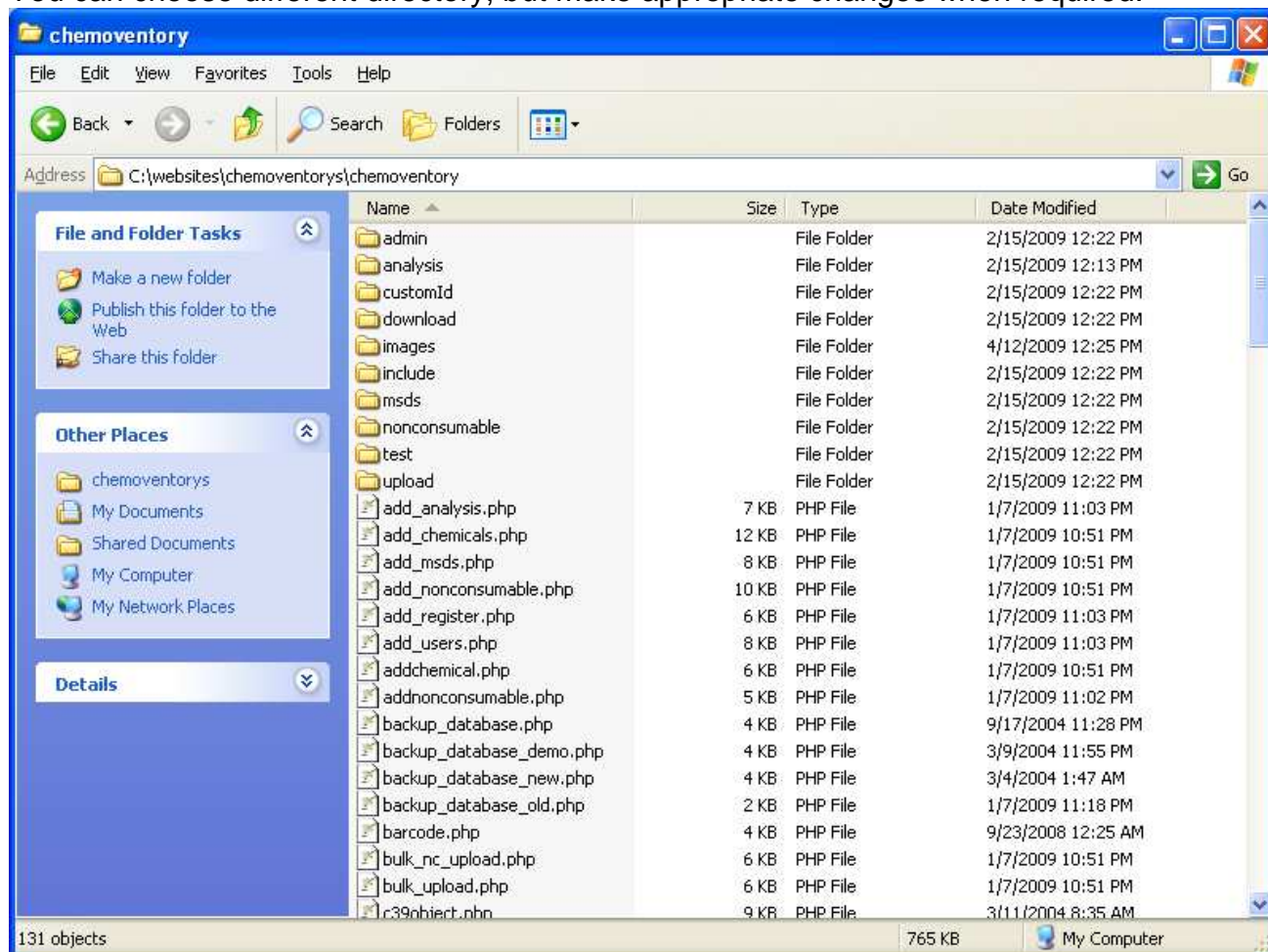


II. APACHE 2.2 Installation in windows XP


1. Open the package




2. Accept conditions

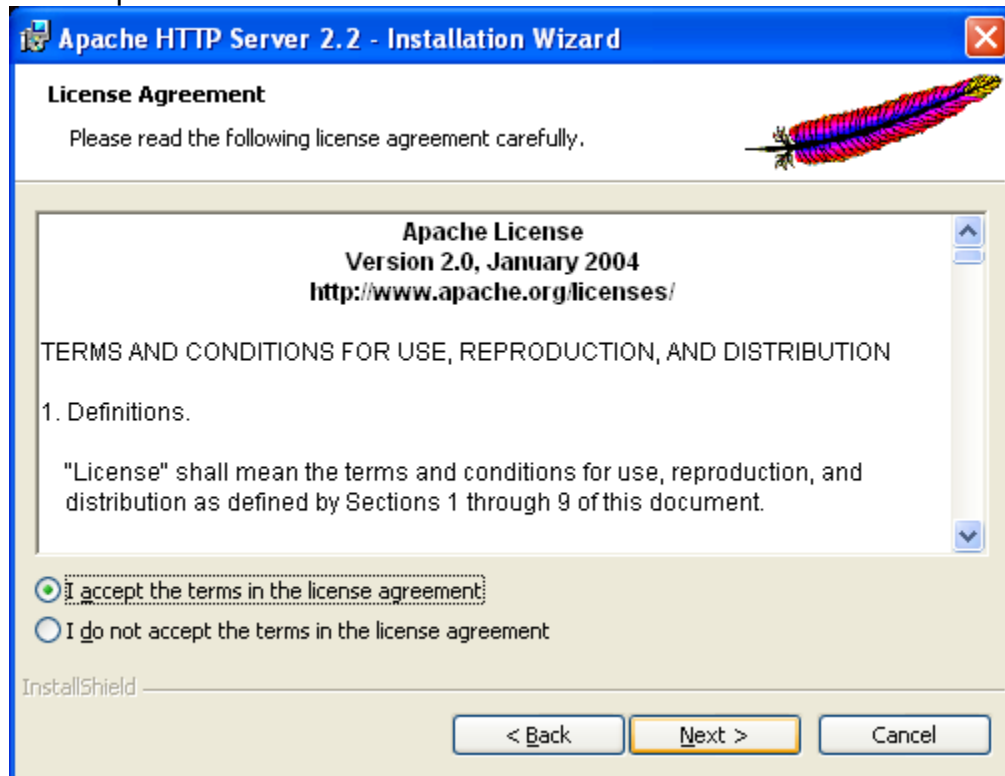


3. Enter your host details

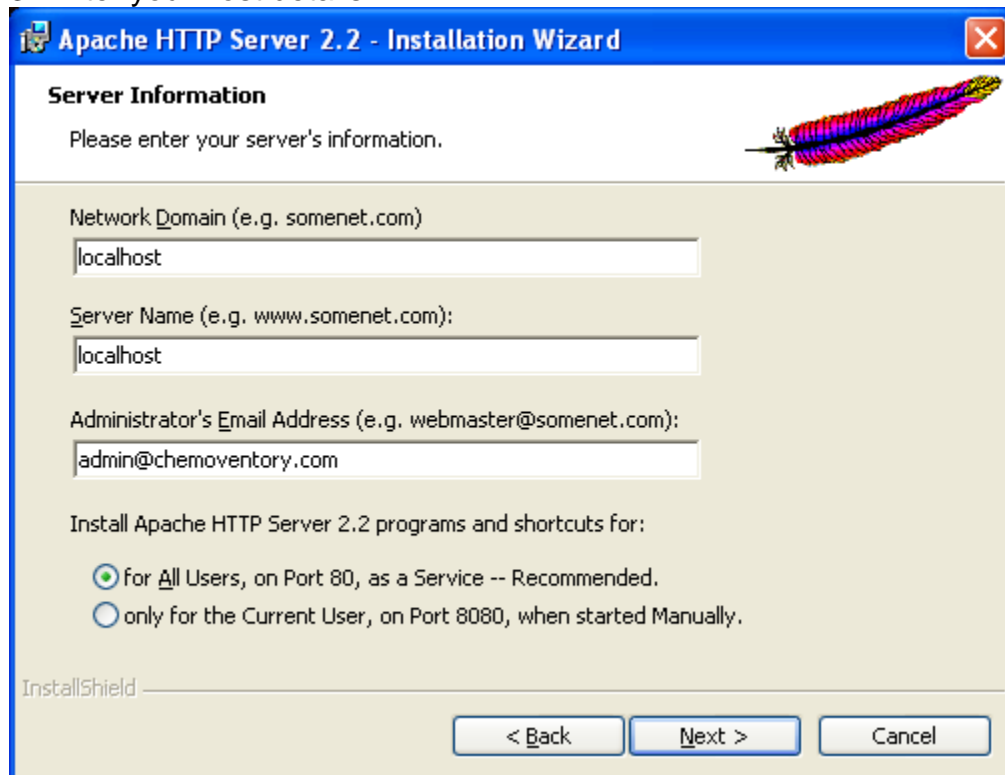


4. Choose the default

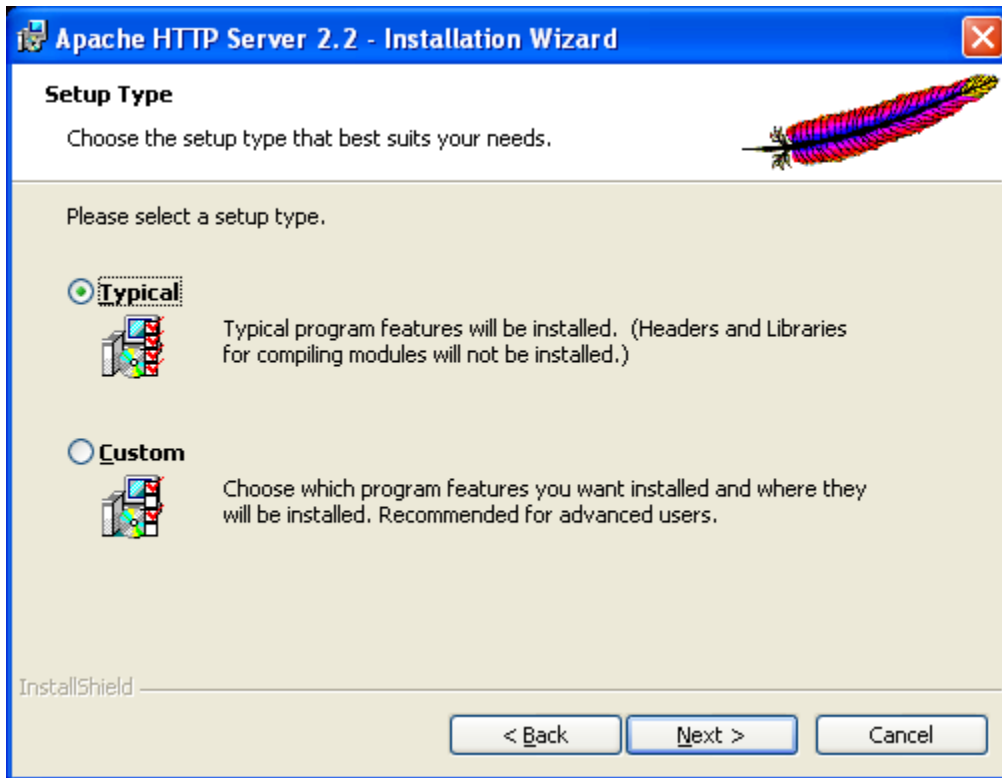


5. Choose the default destination folders

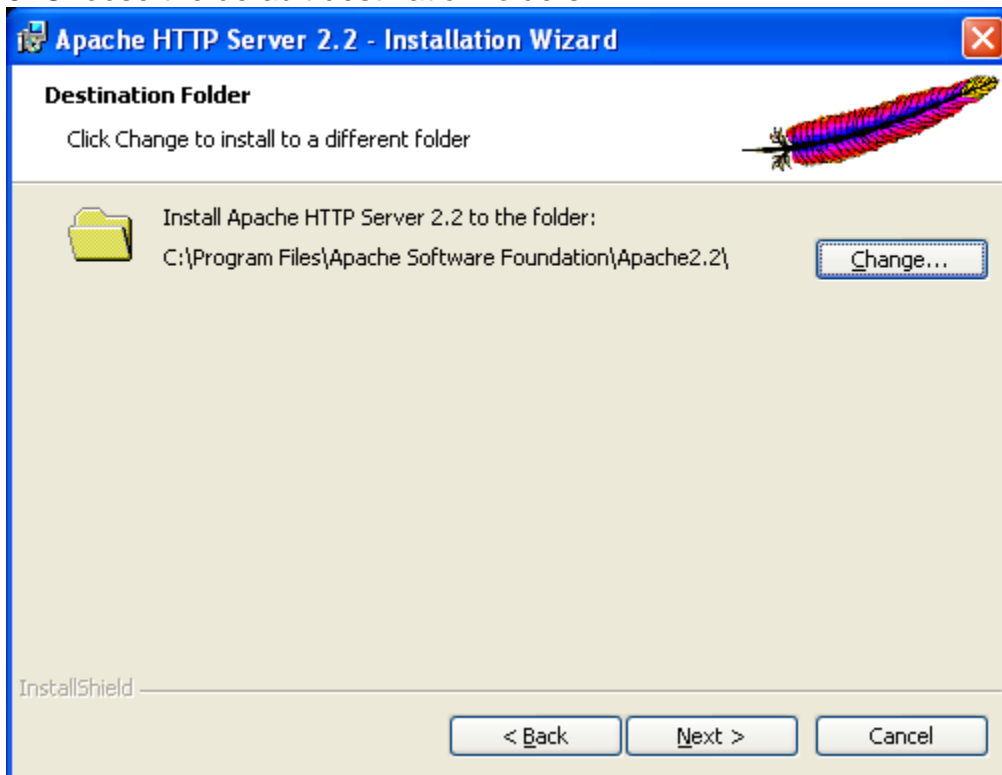

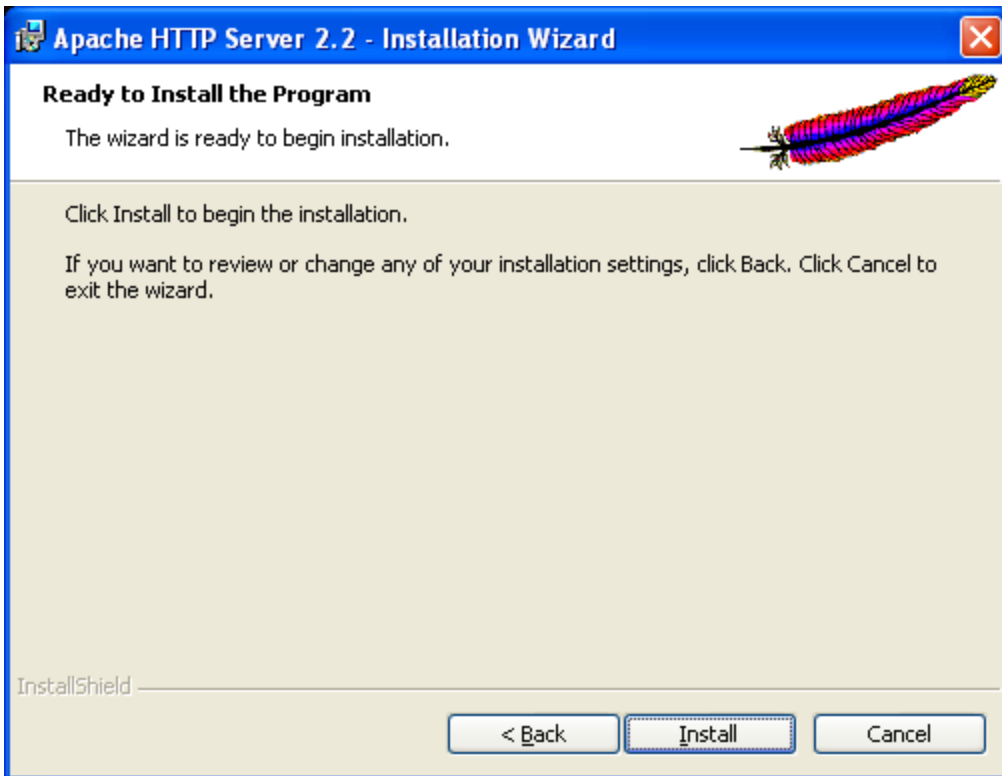


6. You should see the message "Apache2.2 service successfully installed"
If you see errors, you should remove and reinstall apache.

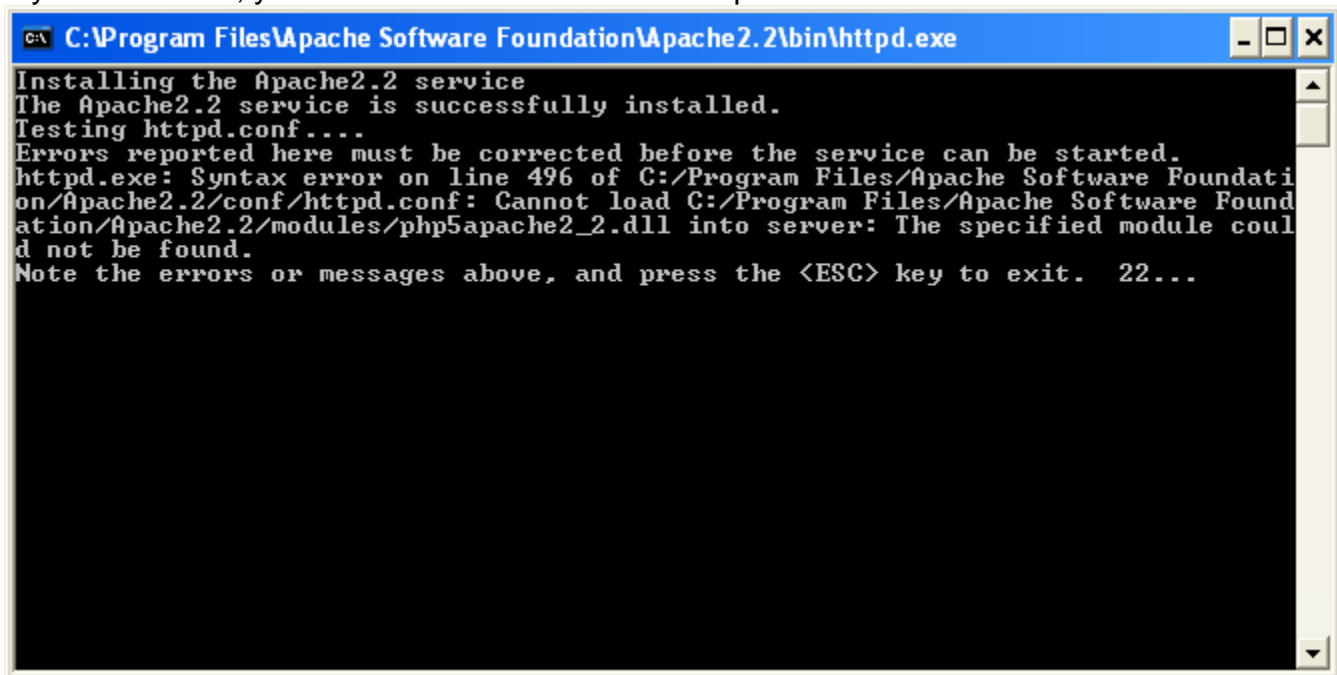

Note: The above error is due to prior apache installation.

Test your installation, before proceeding to next step. In browser visit <http://localhost/>
You should see something in your web page (sometime as shown below or Apache document page)

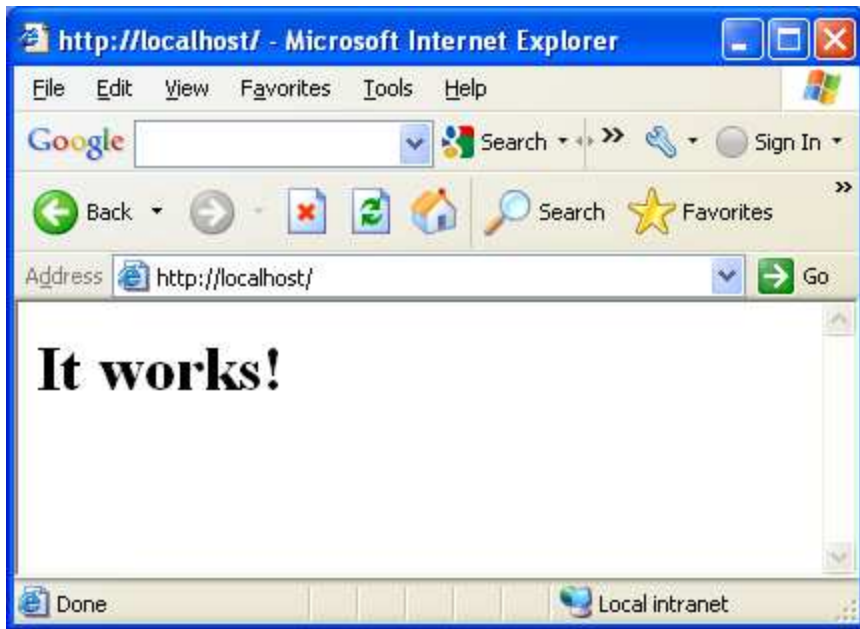


III. PHP 5.2 Installation


Download PHP Binary with installer. <http://www.php.net/downloads.php>




2. Choose the custom folder, instead of default folder

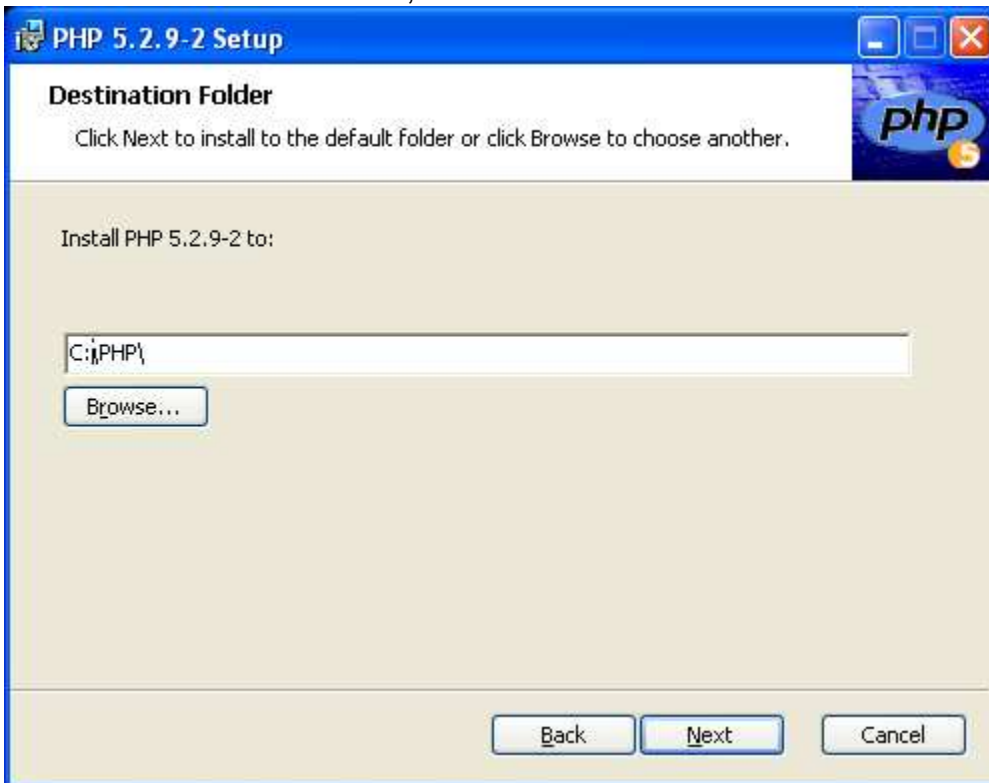



3. Choose your web server (apache2.2)


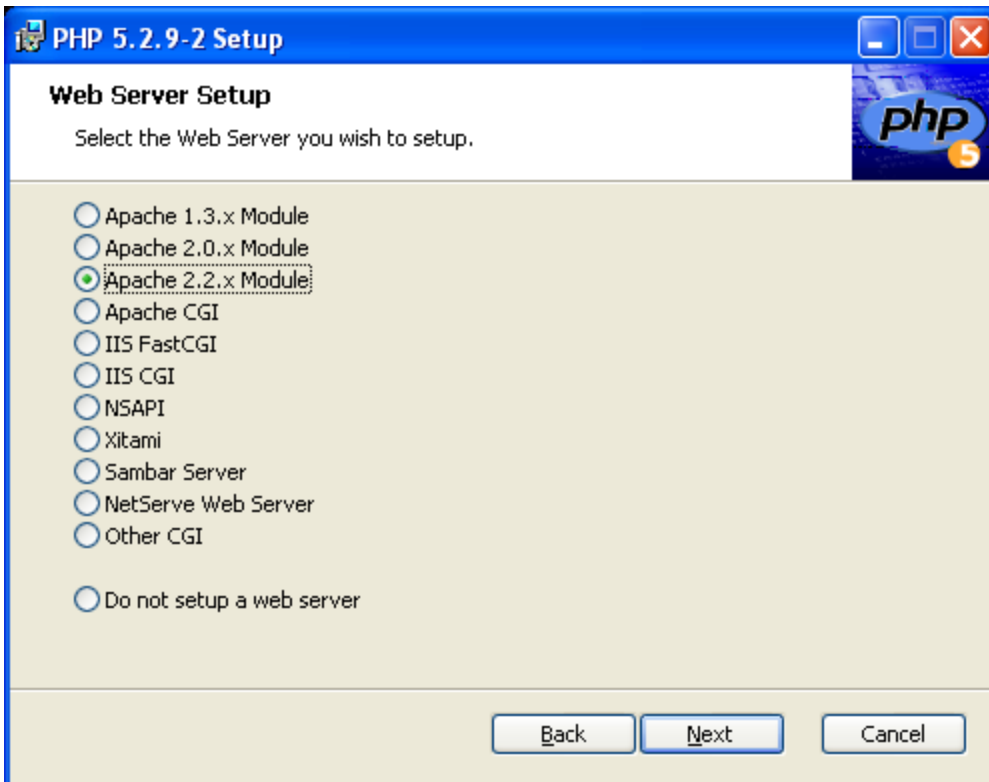


4. You should specify the apache configuration file location httpd.conf. It usually resides in conf directory.

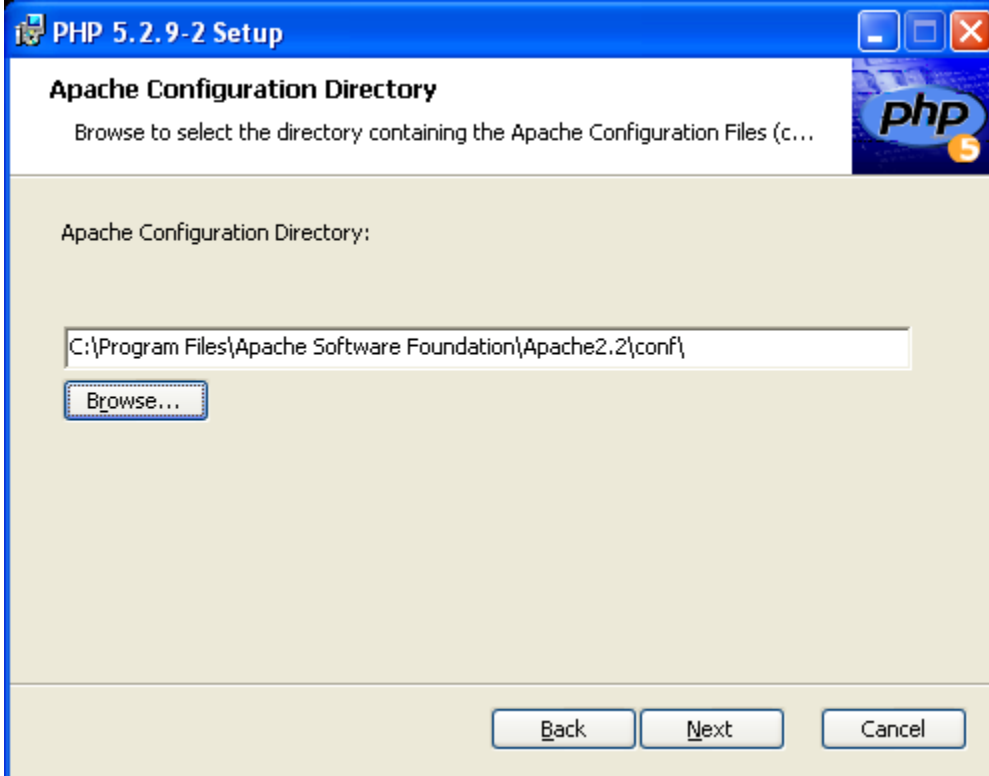

5. You should choose the following three extensions for chemoventory:
GD2, MySQL and PEAR (PEAR may not be required in many cases depending on email system).

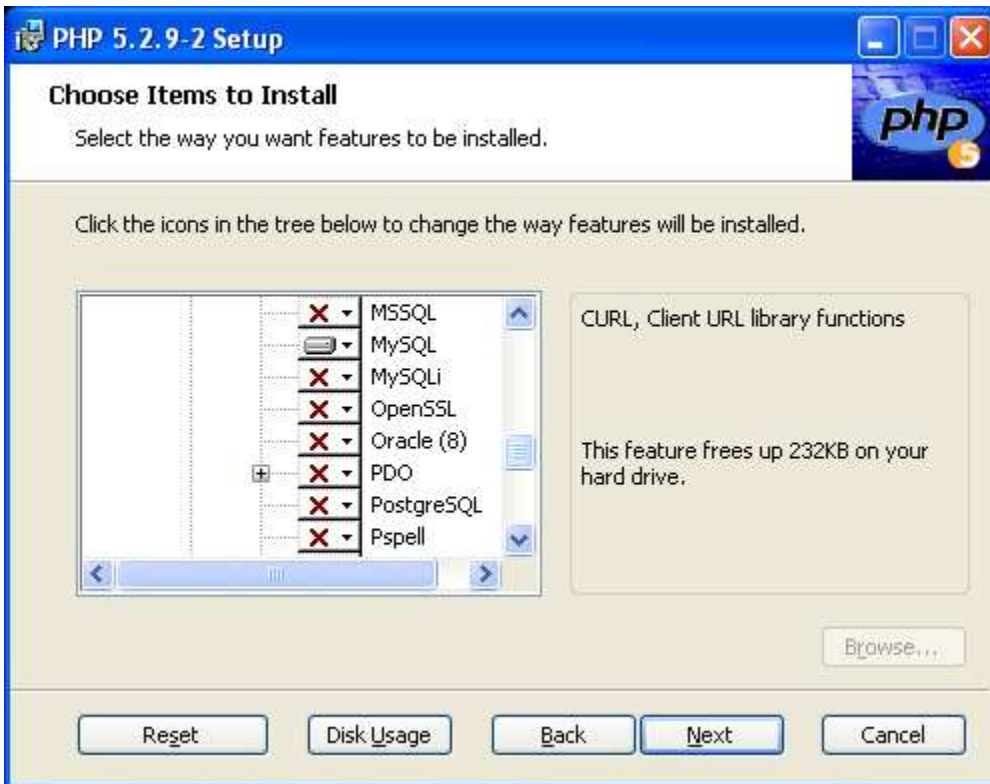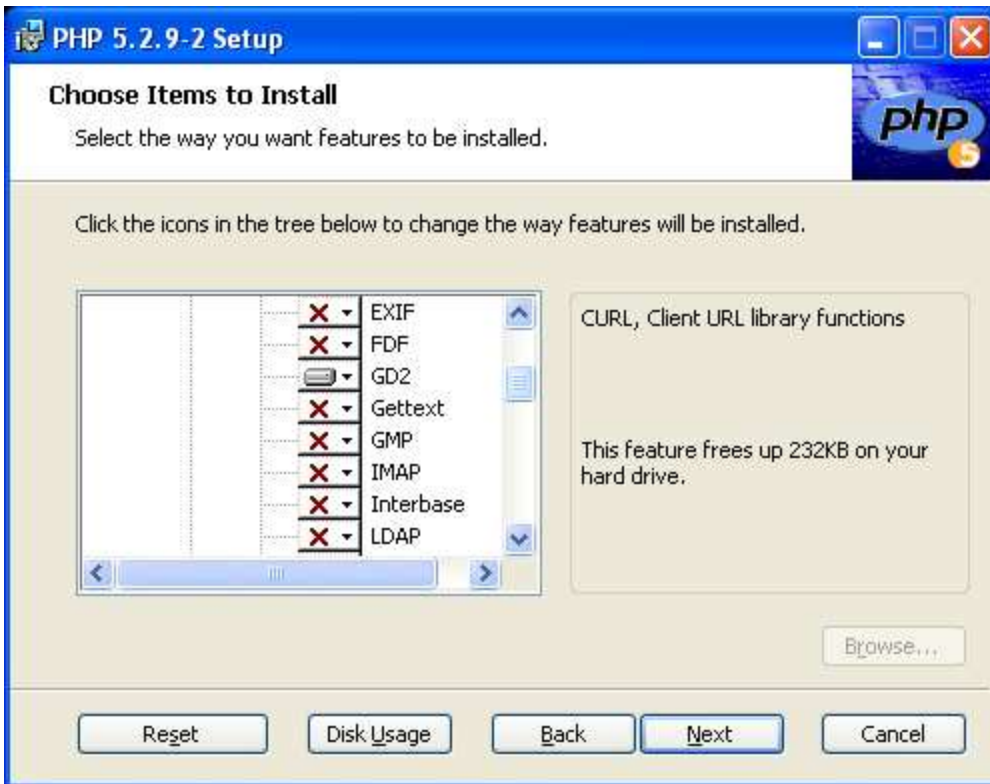

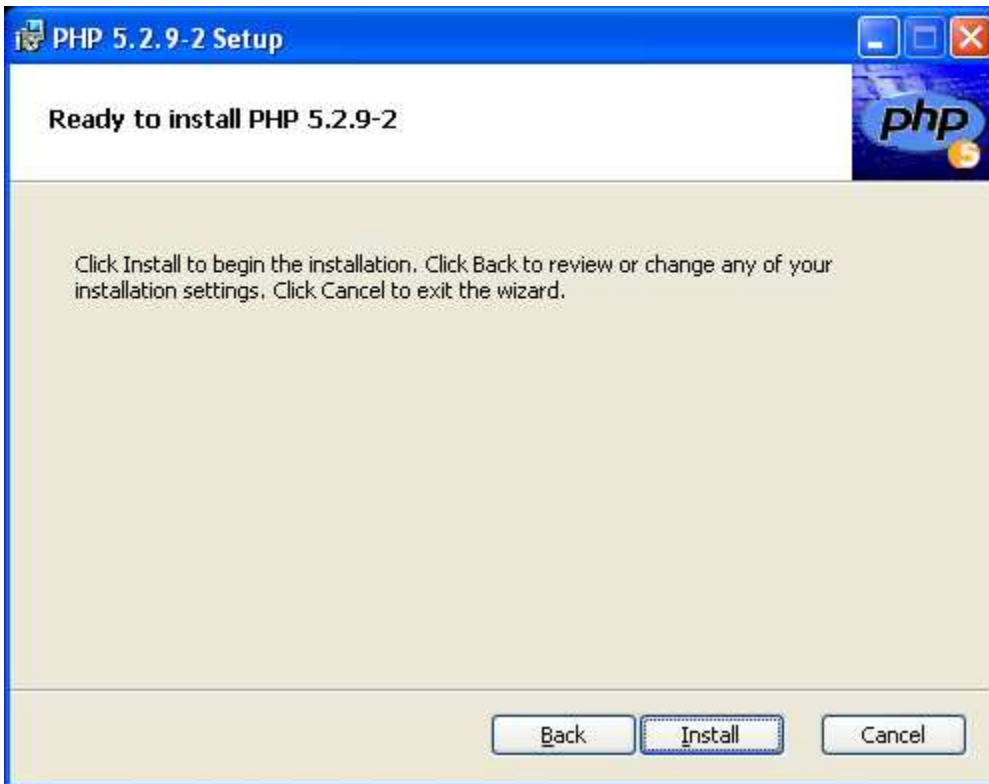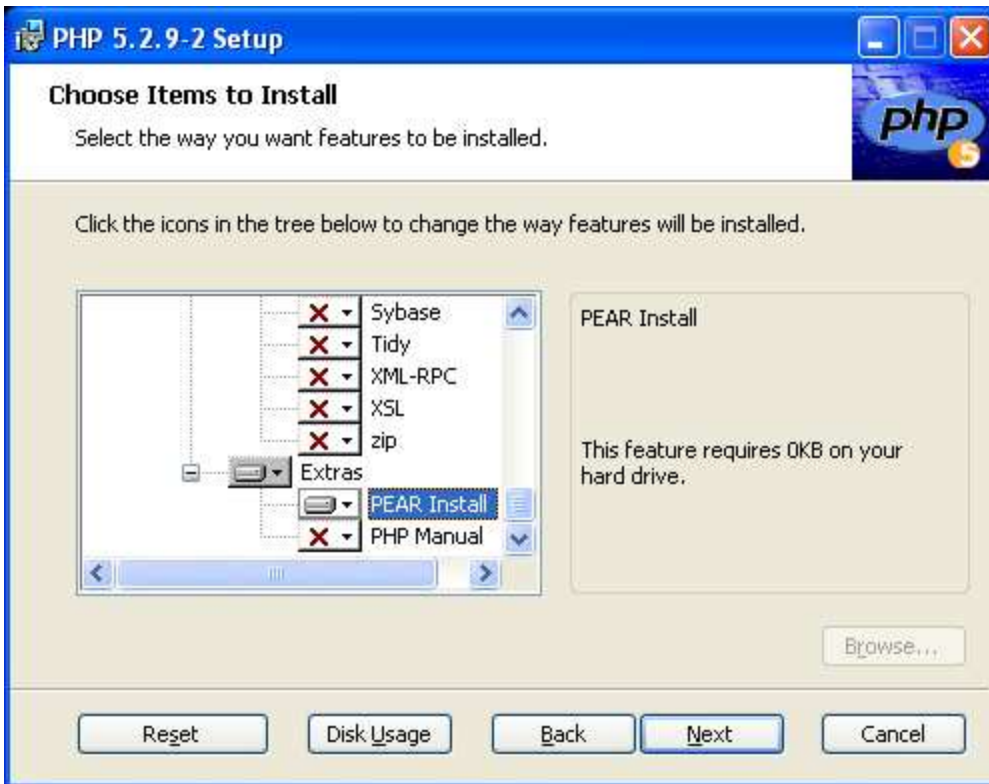

6. Complete the installation



7. Additional configurations in apache's httpd.conf file.

Add the bold one to have php as one of the default file.

```
<IfModule dir_module>  
    DirectoryIndex index.html index.php  
</IfModule>
```

8. Add the entire one to include chemoventory directory as

```
Alias /chemoventory/ "C:/chemoventory/"  
<Directory "C:/chemoventory">  
Options Indexes FollowSymlinks MultiViews  
DirectoryIndex index.html index.html.var index.php  
AllowOverride Indexes  
Order allow,deny  
Allow from all  
</Directory>
```

9. At the end of httpd.conf make sure you have these lines:

```
AddType application/x-httpd-php .php  
#BEGIN PHP INSTALLER EDITS - REMOVE ONLY ON UNINSTALL  
PHPIniDir "C:/PHP/"  
LoadModule php5_module "C:/PHP/php5apache2_2.dll"  
#END PHP INSTALLER EDITS - REMOVE ONLY ON UNINSTALL
```


10. Test PHP installation, by visiting the <http://localhost/chemoventory/test/test.php>

phpinfo() - Microsoft Internet Explorer

File Edit View Favorites Tools Help Google Search Sign In

Address http://localhost/websites/chemoventorys/chemoventory/test/test.php Go



PHP Version 5.2.9-2





System	Windows NT CHEMO 5.1 build 2600
Build Date	Apr 9 2009 08:22:37
Configure Command	<pre>cscript/nologo configure.js "--enable-snapshot-build" "--enable-debug-pack" "--with-snapshot-template=d:\php-sdk\snap_5_2\vc6x86\template" "--with-php-build=d:\php-sdk\snap_5_2\vc6x86\php_build" "--with-pdo-oci=D:\php-sdk\oracle\instantclient10\sdk_shared" "--with-oci8=D:\php-sdk\oracle\instantclient10\sdk_shared"</pre>
Server API	Apache 2.0 Handler
Virtual Directory Support	enabled
Configuration File (php.ini) Path	D:\WINDOWS
Loaded Configuration File	C:\php\php.ini
Scan this dir for additional .ini files	(none)
additional .ini files parsed	(none)
PHP API	20041225
PHP Extension	20060613
Zend Extension	220060519
Debug Build	no
Thread Safety	enabled
Zend Memory Manager	enabled



Done Local intranet

IV. Zend Encoder

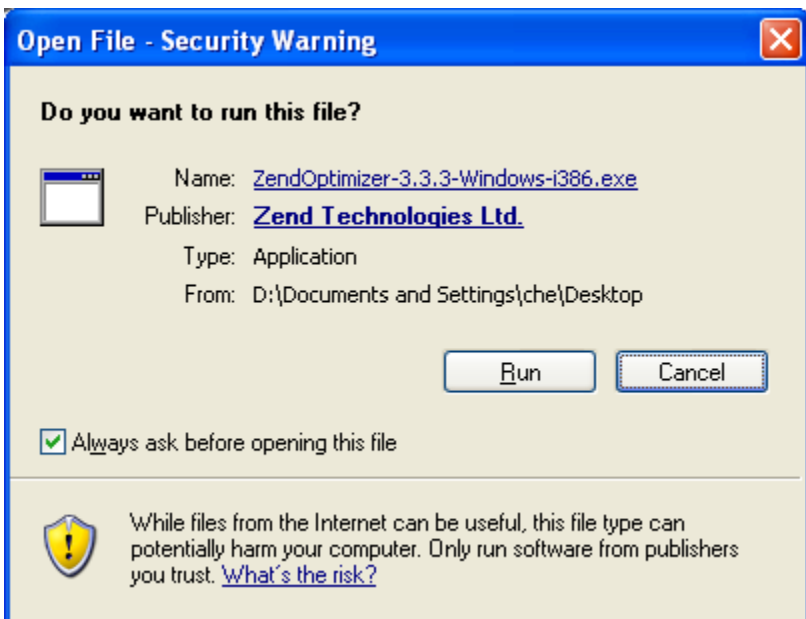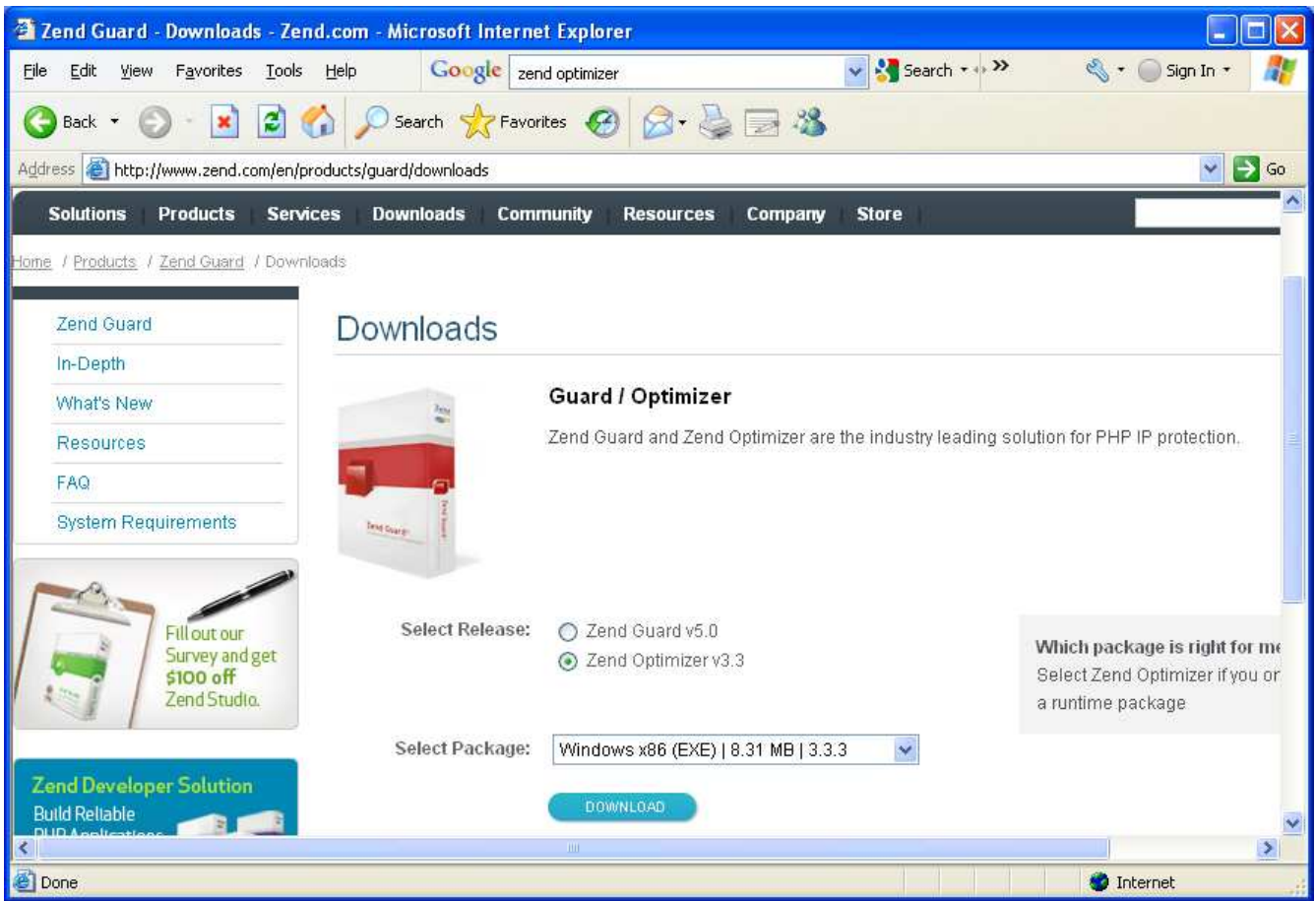

Select option of apache2.2

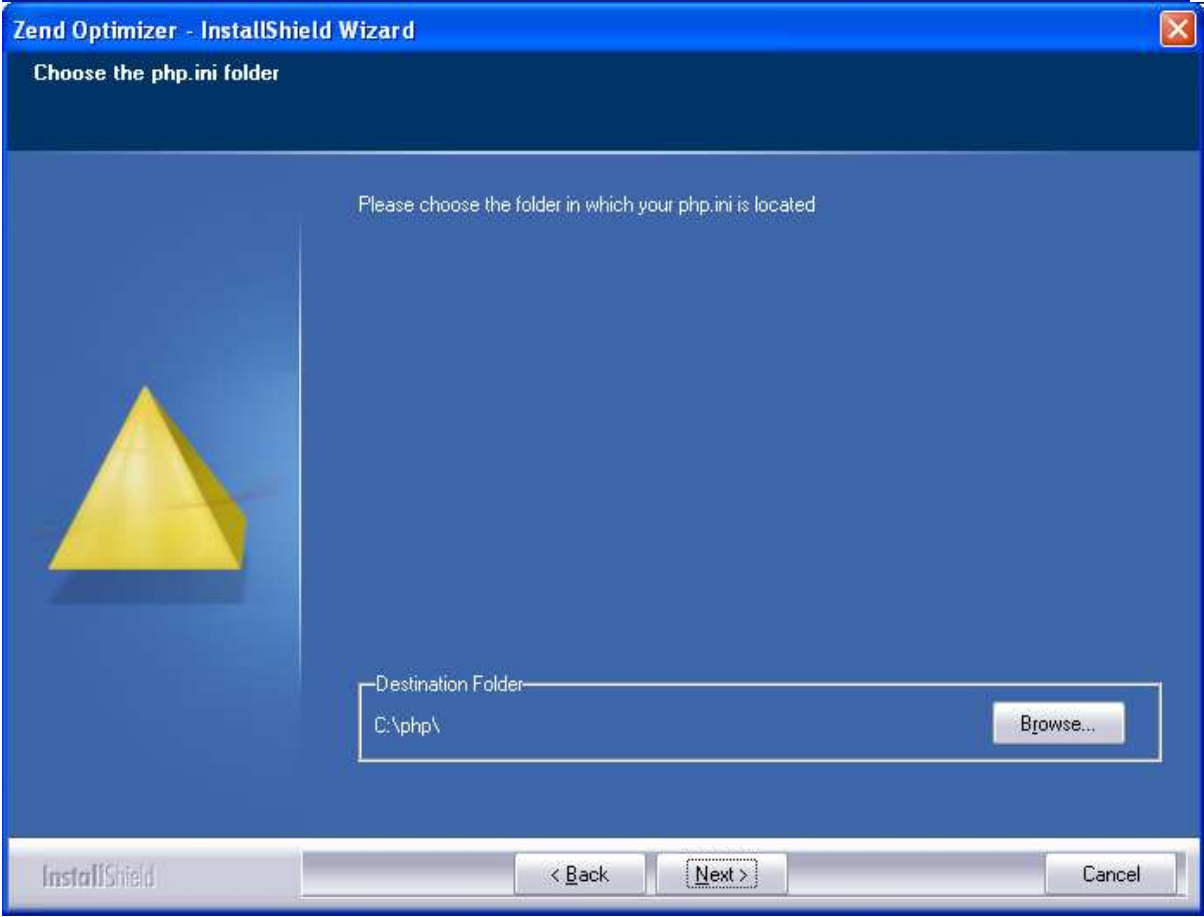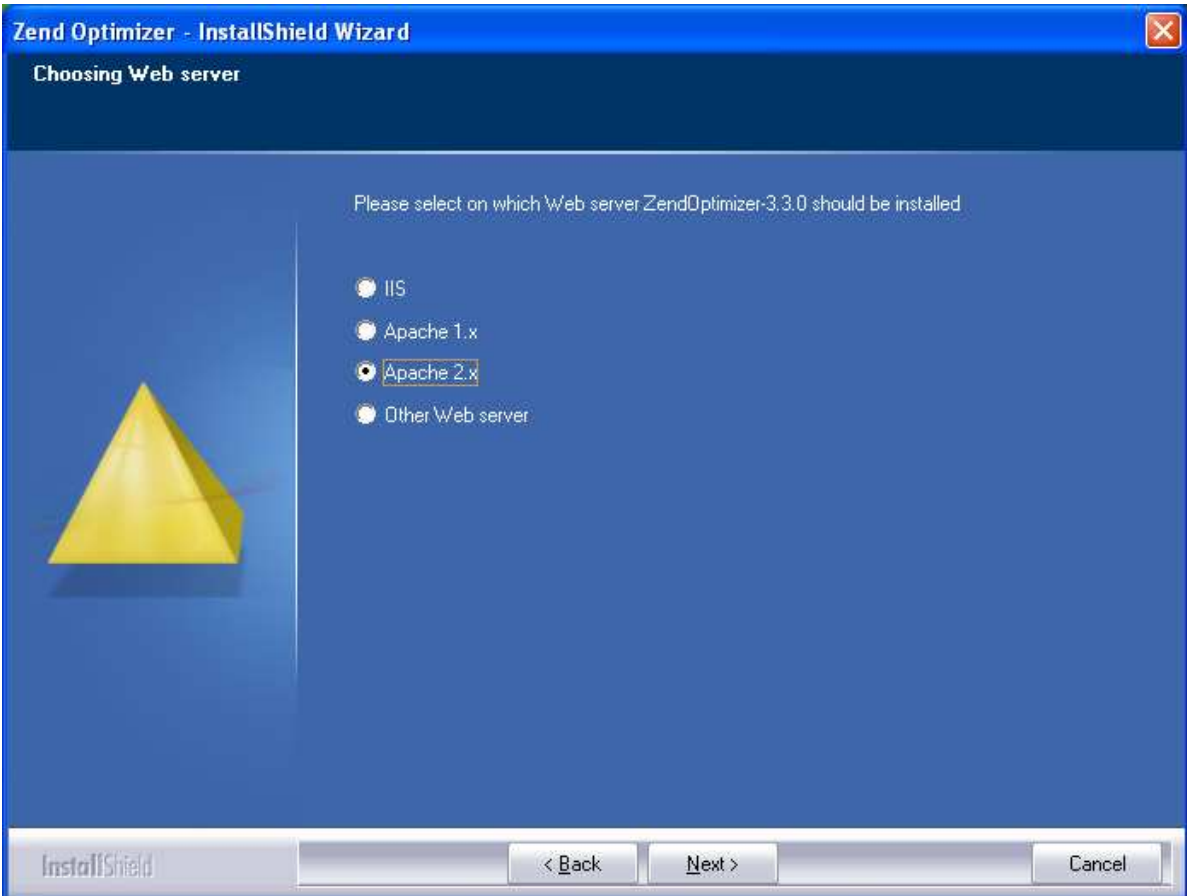

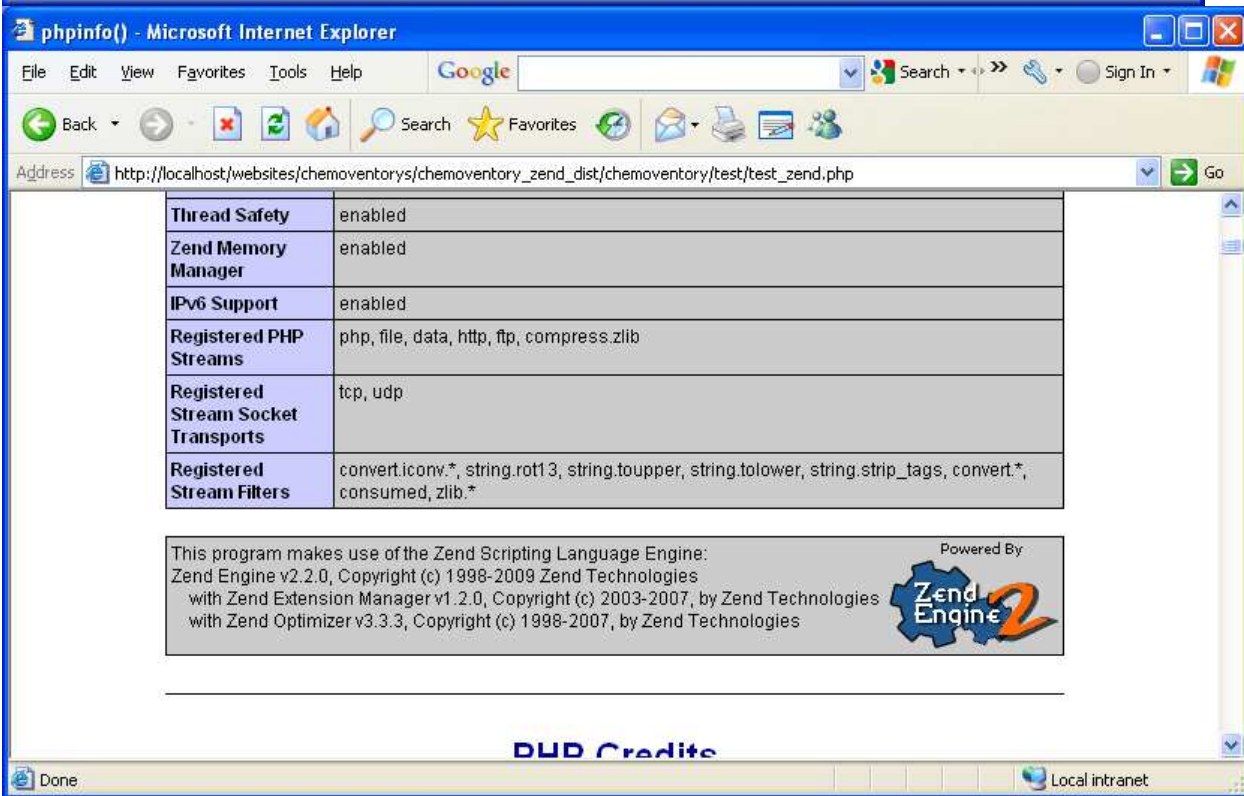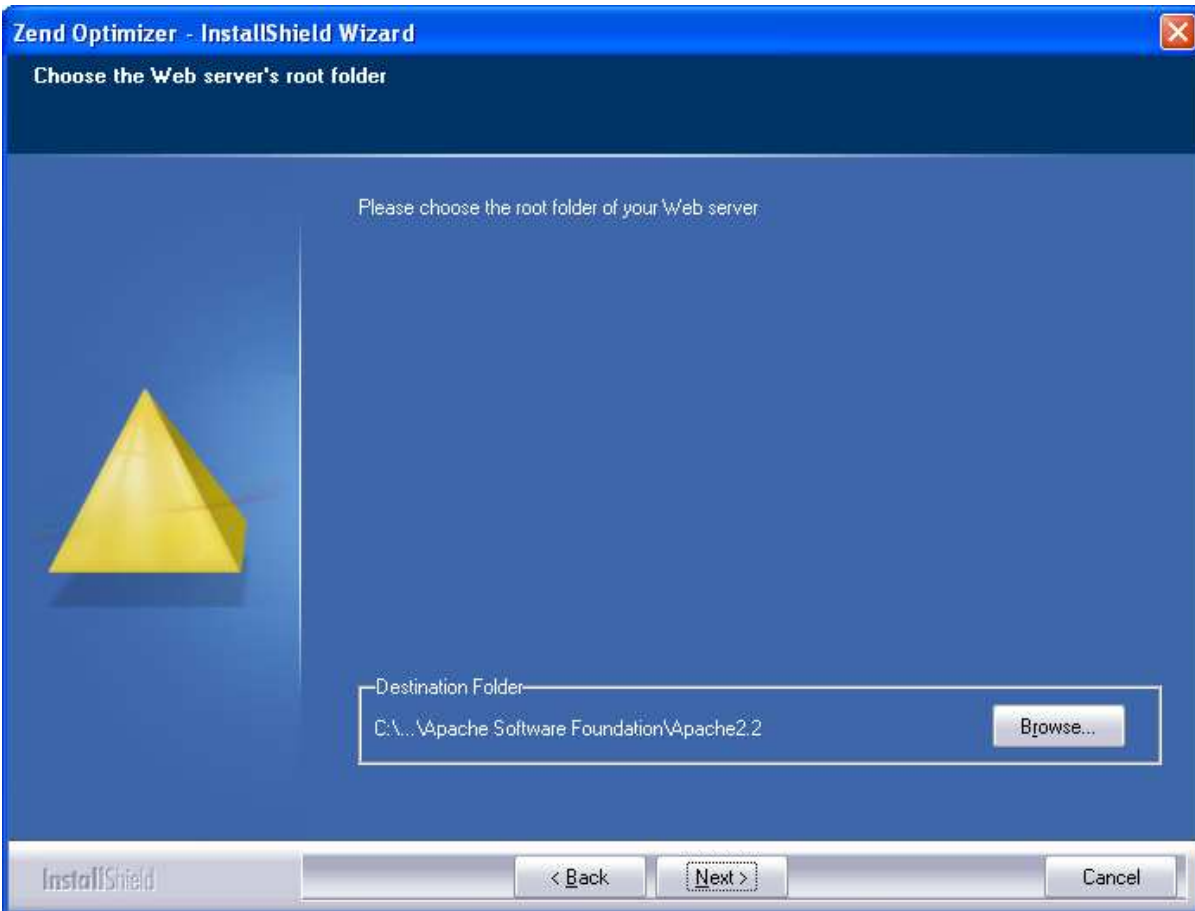



At the end of php.ini (c:/php) you should see

```
[Zend]
zend_extension_manager.optimizer_ts="D:\Program Files\Zend\ZendOptimizer-3.3.0\lib\Optimizer-3.3.0"
zend_extension_ts="D:\Program Files\Zend\ZendOptimizer-3.3.0\lib\ZendExtensionManager.dll"
```




V. MySQL Installation


Here are the steps to install MySQL




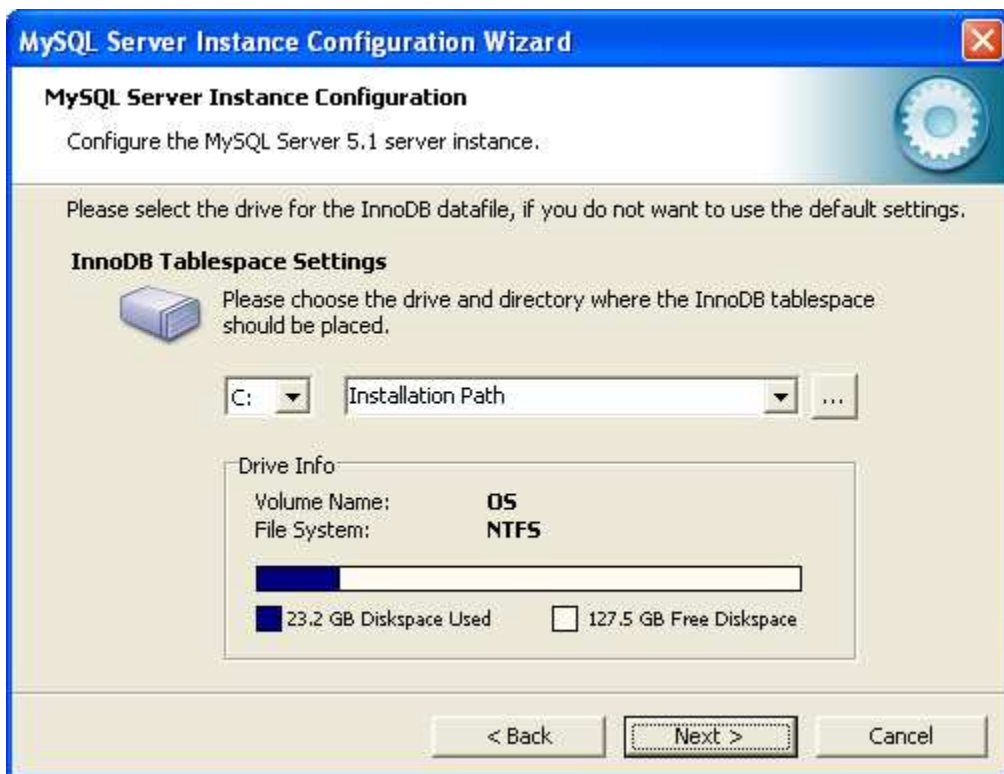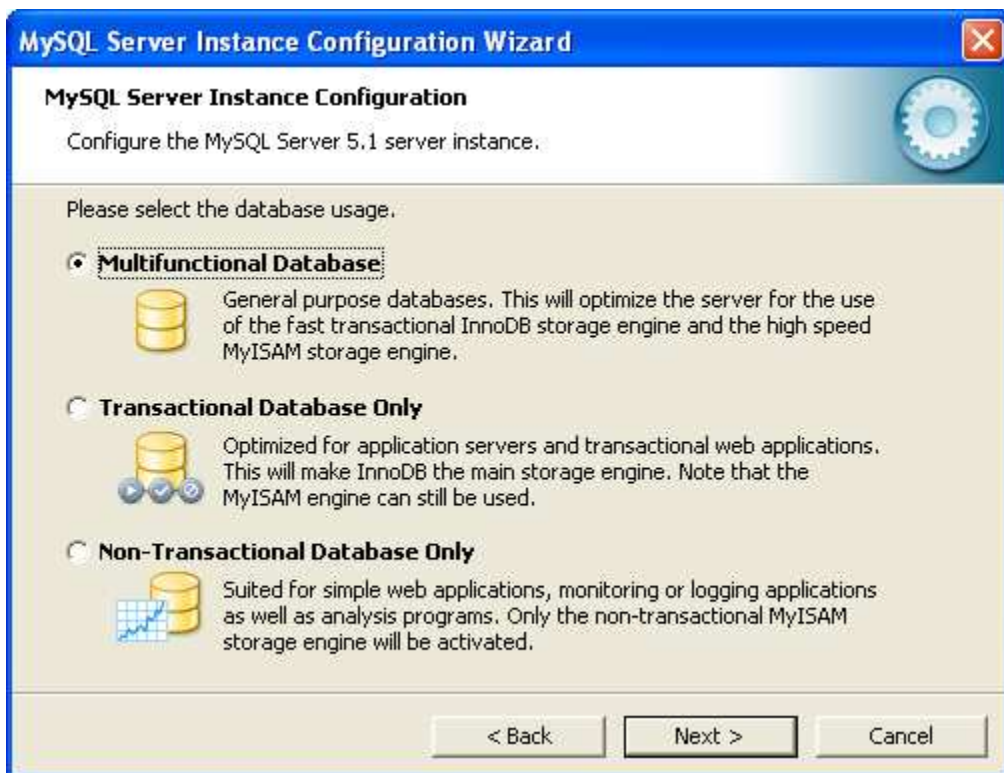

Remove the strict mode option

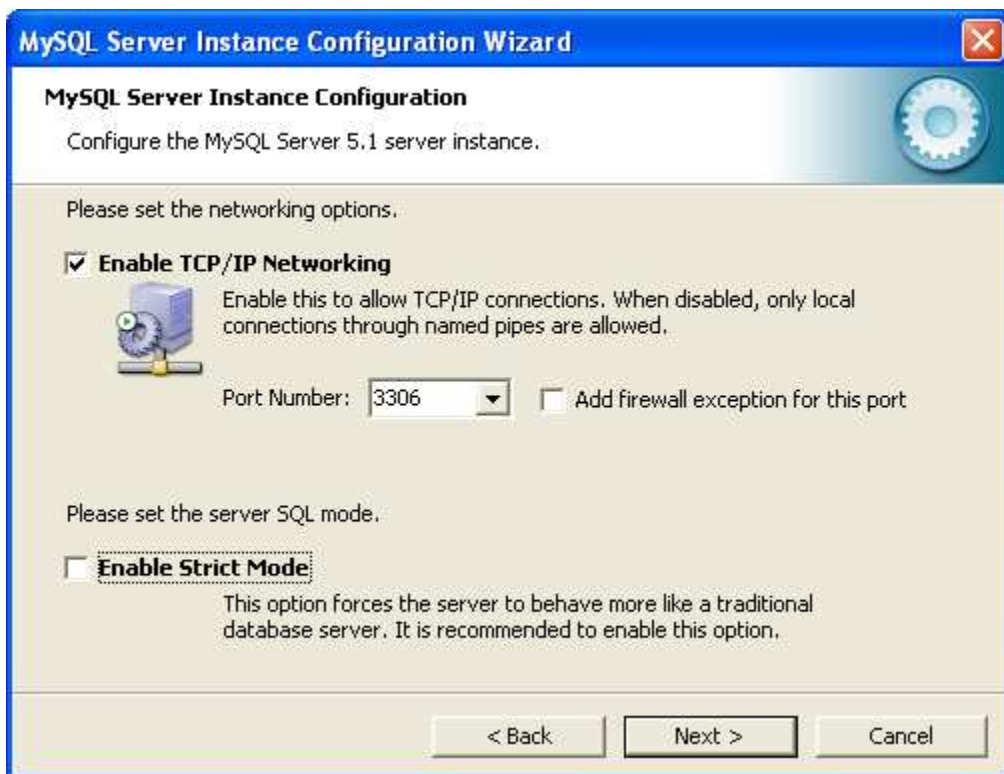
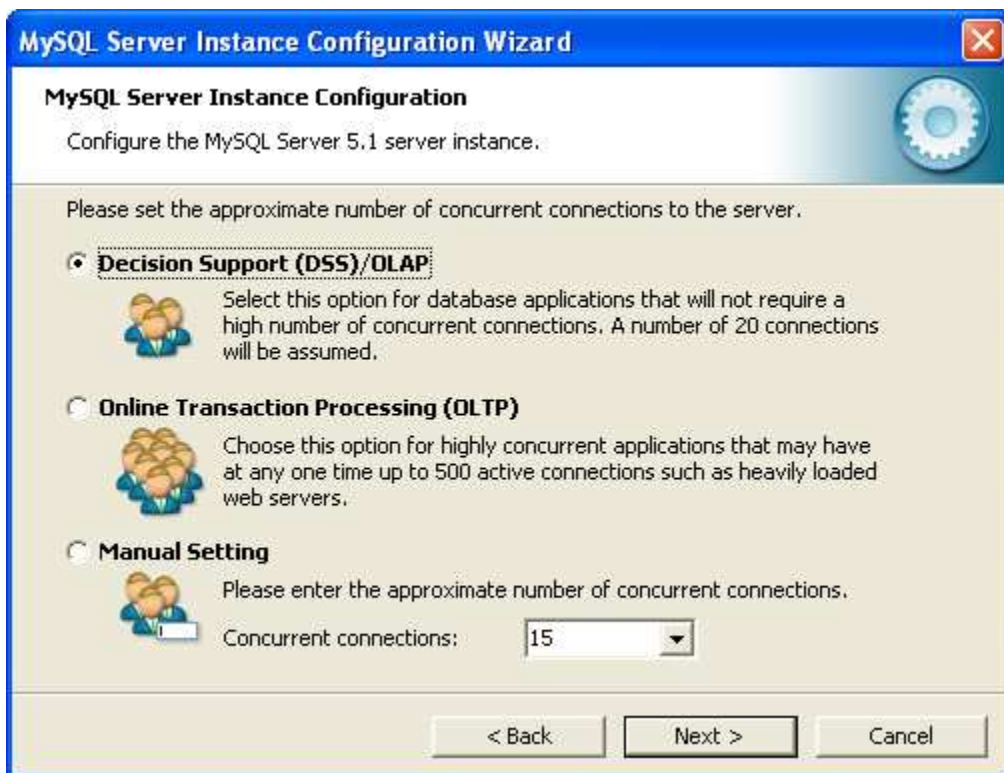


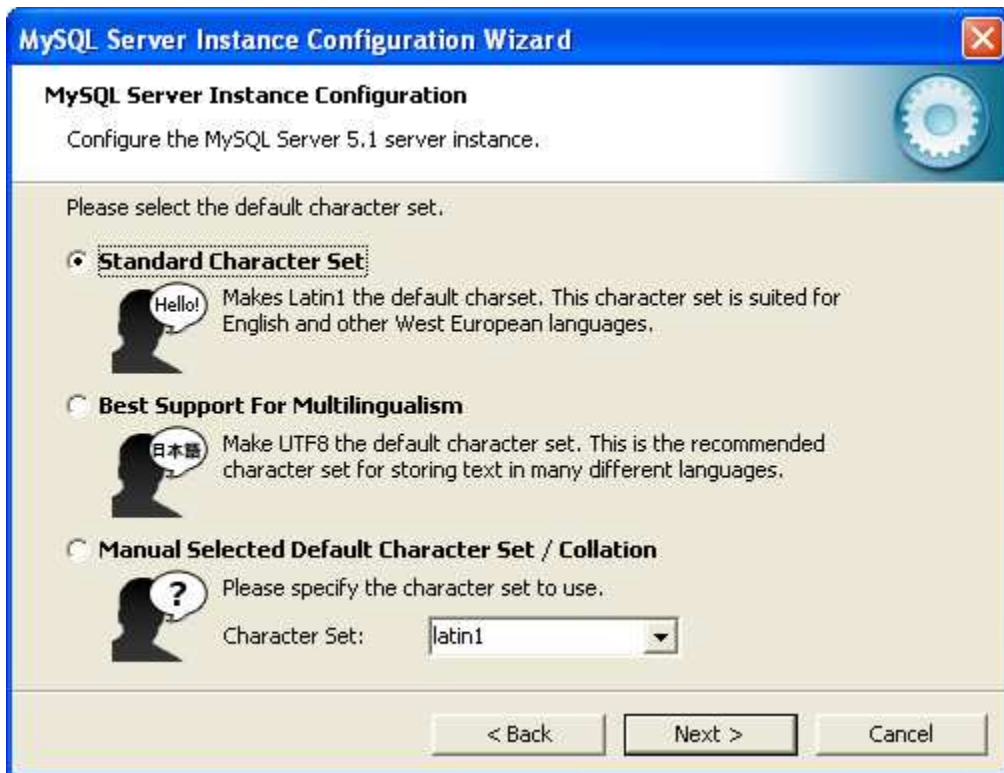

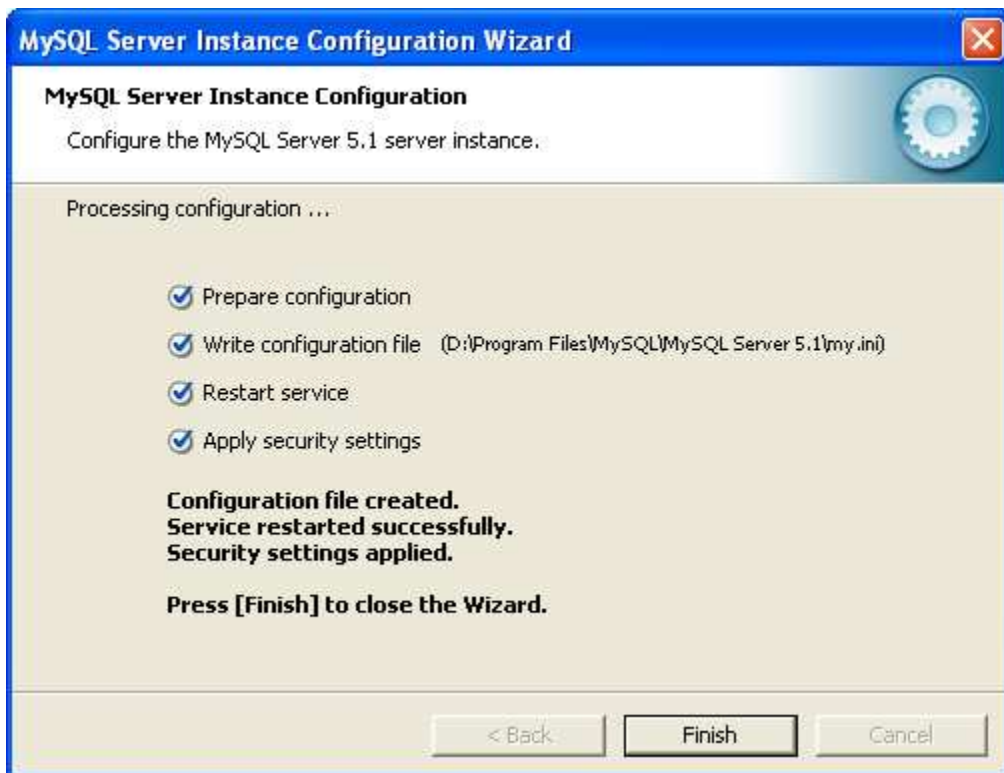

Install GUI interface for MySQL

Download Heidi MySQL from <http://www.heidisql.com/download.php> (with Installer)

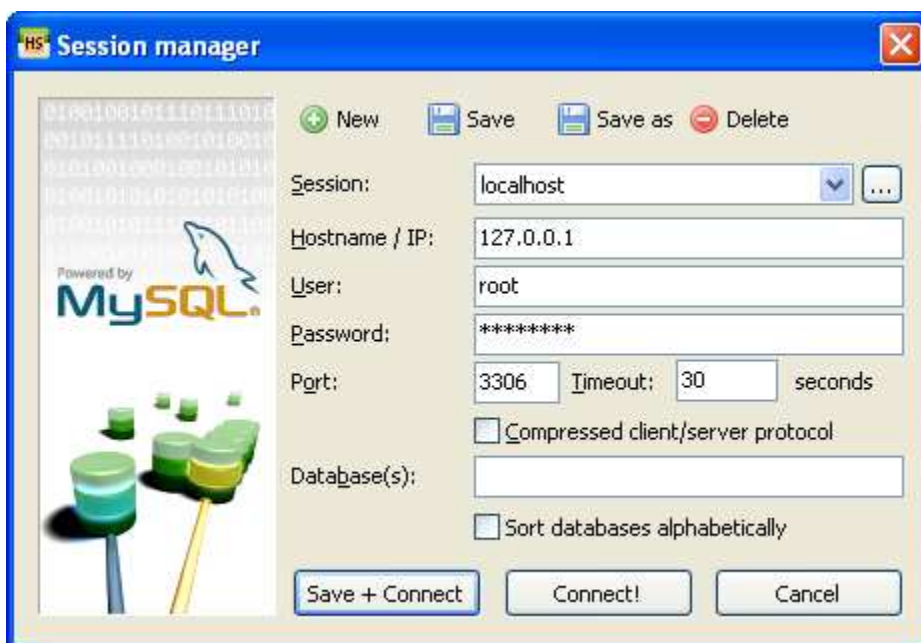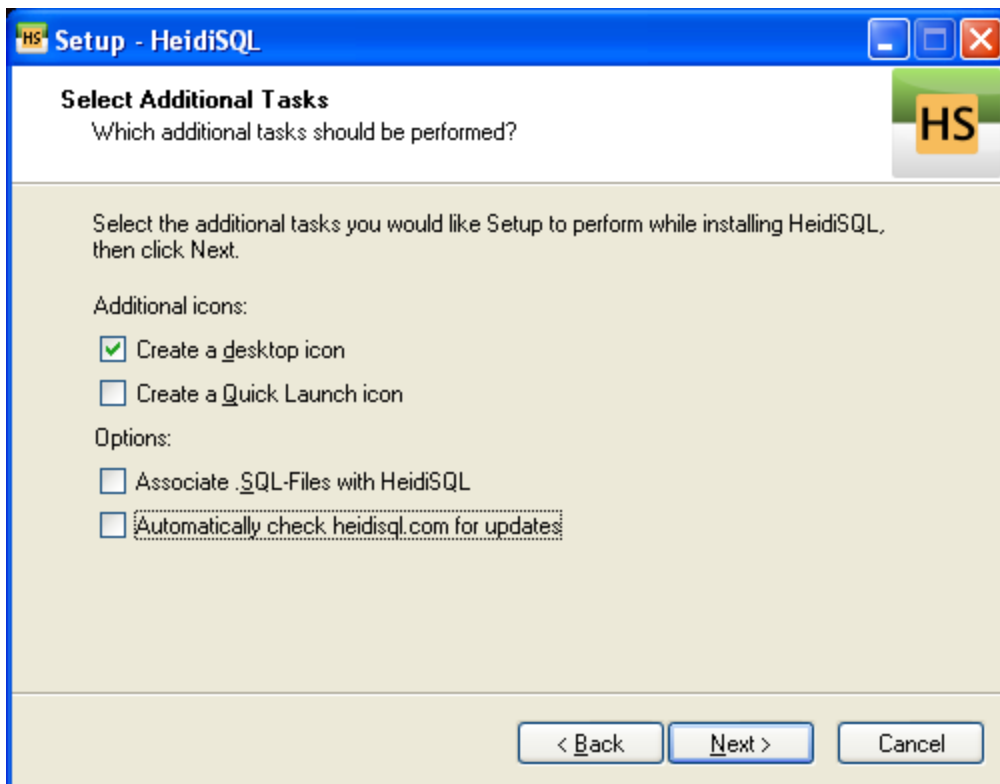

localhost - /chemoventory

File Edit Tools Import Export Window Help

Host: 127.0.0.1 Database: chemoventory Query

Name	Rows	Size	Created	Updated	Engine	Comment	Type
chemo_analysis	0	1.0 KB	2009-05-15 23:42:34	2009-05-15 23:42:34	MyISAM		BASE...
chemo_audit	0	32.0 KB	2009-05-15 23:42:34		InnoDB		BASE...
chemo_buildings	3	3.1 KB	2009-05-15 23:42:34	2009-05-15 23:43:04	MyISAM		BASE...
chemo_category	15	4.3 KB	2009-05-15 23:42:34	2009-05-15 23:43:04	MyISAM		BASE...
chemo_chemical	1	64.0 KB	2009-05-15 23:42:34		InnoDB		BASE...
chemo_column	6	1.1 KB	2009-05-15 23:42:34	2009-05-15 23:43:04	MyISAM		BASE...
chemo_configur...	3	3.1 KB	2009-05-15 23:42:34	2009-05-15 23:43:05	MyISAM		BASE...
chemo_containers	49	5.5 KB	2009-05-15 23:42:34	2009-05-15 23:44:11	MyISAM		BASE...
chemo_dailyjob	1	2.0 KB	2009-05-15 23:42:34	2009-05-15 23:43:30	MyISAM		BASE...
chemo_deleted	0	1.0 KB	2009-05-15 23:42:34	2009-05-15 23:42:34	MyISAM		BASE...
chemo_divisions	6	4.1 KB	2009-05-15 23:42:34	2009-05-15 23:43:05	MyISAM		BASE...
chemo_functions	49	11.0 KB	2009-05-15 23:42:34	2009-05-15 23:44:11	MyISAM		BASE...
chemo_locations	0	1.0 KB	2009-05-15 23:42:34	2009-05-15 23:42:34	MyISAM		BASE...
chemo_logactivity	0	32.0 KB	2009-05-15 23:42:34		InnoDB		BASE...
chemo_logs	2	16.0 KB	2009-05-15 23:42:34		InnoDB		BASE...
chemo_manufa...	3	4.1 KB	2009-05-15 23:42:34	2009-05-15 23:43:06	MyISAM		BASE...
chemo_msds	0	2.0 KB	2009-05-15 23:42:34	2009-05-15 23:42:34	MyISAM		BASE...
chemo_nc_cate...	6	4.1 KB	2009-05-15 23:42:34	2009-05-15 23:43:06	MyISAM		BASE...
chemo_nc_dele...	0	1.0 KB	2009-05-15 23:42:34	2009-05-15 23:42:34	MyISAM		BASE...
chemo_noncons...	0	16.0 KB	2009-05-15 23:42:34		InnoDB		BASE...
chemo_owner	0	48.0 KB	2009-05-15 23:42:34		InnoDB		BASE...
chemo_properties	0	2.0 KB	2009-05-15 23:42:34	2009-05-15 23:42:34	MyISAM		BASE...
chemo_recycle	0	1.0 KB	2009-05-15 23:42:34	2009-05-15 23:42:34	MyISAM		BASE...
chemo_reports	0	1.0 KB	2009-05-15 23:42:34	2009-05-15 23:42:34	MyISAM		BASE...
chemo_roles	88	7.4 KB	2009-05-15 23:42:34	2009-05-15 23:44:11	MyISAM		BASE...
chemo_temp_in...	0	16.0 KB	2009-05-15 23:42:34		InnoDB		BASE...
chemo_temp_nc	0	16.0 KB	2009-05-15 23:42:35		InnoDB		BASE...

```

7 USE `chemoventory` ;
8 SHOW TABLES FROM `information_schema` ;
9 (SELECT TABLE_NAME AS `Name`, TABLE_TYPE AS `Type`, ENGINE AS `Engine`, VERSION AS `Version`, ROW_FORMAT AS `Row_format`, TABLE_ROWS AS `Rows`, AVG_ROW_LENGTH AS `Avg_row_leng

```

chemoventory: 33 table(s) Connected: 00:00:13 MySQL 5.1.34 Uptime: 0 days, 04:09:29 Ready.